



MEDIGENIC®

HYGIENIC KEYBOARDS AND MICE
WITH A HYGIENIC WARNING SYSTEM, DISINFECTION IN SECONDS
SAME TYPING EXPERIENCE AS A NORMAL KEYBOARD

A close-up photograph of a person wearing green scrubs and a blue surgical cap, washing their hands under a running faucet in a stainless steel sink. The water is clear and flowing, splashing slightly. The background is blurred, showing a typical clinical setting.

SAFE & CLEAN

Keyboards and mice have become vital in medical environments, operating rooms and at nursing stations. However, traditional input devices are usually unsuitable for medical use due to hygiene reasons. The MEDIGENIC keyboards and mice are simple and safe to clean and make it easier to comply with the required hygiene standards without foregoing ease of use.

MEDIGENIC KEYBOARDS & MICE

The MEDIGENIC keyboards and mice facilitate compliance with specified hygiene standards without foregoing ease of use. The specially developed design provides an excellent typing experience with tangible button movement. User operating speeds remain at a high level and there is no need to fear comparison with traditional keyboards. MEDIGENIC keyboards and mice are a reliable choice in all environments where hygiene and health are top priorities.



Hygienically clean in seconds

The flat, ergonomic keyboard design enables fast cleaning with all hospital-approved disinfectants. The moisture-resistant surface based on silicone is completely impermeable to liquids, latex-free and extremely durable.

Optimal typing feel

The Medigenic keyboards provide solid tactile sensory feedback when operated. This enables fast and reliable data input. The 3D surface of the key panel simulates the haptics of a keyboard with standard buttons and thereby makes it easier to position the fingers.



MEDIGENIC FAMILY OVERVIEW




MEDIGENIC ESSENTIAL

Medigenic Essential is an easy-to-clean, cost-effective keyboard ideal for use in medical environments. The special product design makes it easy to clean and disinfect. The hypoallergenic button panel is latex-free and does not contain other irritating materials.



MEDIGENIC COMPLIANCE

Medigenic Compliance offers all of the same benefits as the Essential keyboard and also features a three-level backlight. An integrated cleaning alarm reminds operating staff when cleaning is required at predefined intervals. As a result, the keyboard automatically detects when cleaning is completed.

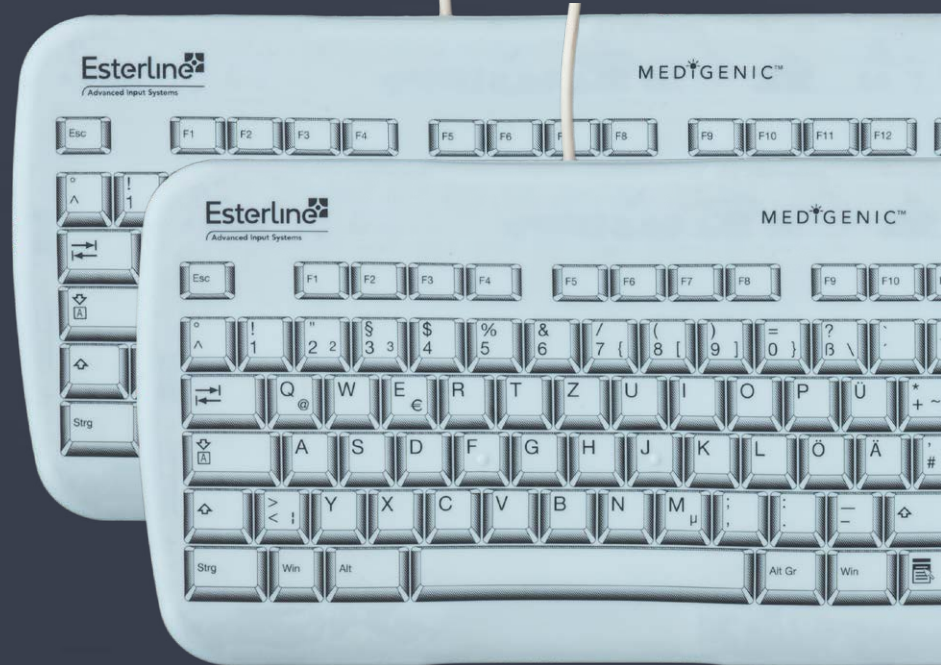

MEDIGENIC COMPACT

Similar to the Compliance keyboard, the Compact keyboard also features a cleaning alarm and a backlight, but has much smaller space requirements. As a result, it is ideally suited for restricted input situations, such as on mobile rounds trolleys.



MEDIGENIC Scroll Mouse

The MEDIGENIC Scroll Mouse is impervious to anti-bacterial sprays and liquid cleaning products due to the fact that it is manufactured using a flexible silicone material. The smooth, latex-free surface is hypoallergenic and does not feature any corners, edges or niches where bacteria can hide. The ergonomic mouse features three touch buttons (left, right, center) and a scroll function.

PERFECT HYGIENE

Patented cleaning alarm

The Compliance and Compact keyboards from Medigenic feature a cleaning alarm that users can individually set. After a defined interval, the keyboard begins to flash, making users aware that cleaning is required. The display disappears automatically when cleaning is completed.

Disinfection without interruption

The keyboard can be temporarily disabled with the push of a button in order to disinfect the device. This prevents incorrect entries. It is not necessary to release the USB connection or shut down the PC. The keyboard is disinfected immediately.

Illuminated touch buttons

The Compliance keyboard is easy and convenient to operate, even in poorly illuminated environments such as those used for ultrasound examinations. This is thanks to integrated button illumination that is adjustable at three levels.

SUMMARY

- Flat, easy-to-clean design - no unhygienic edges and corners
- IP65 protection against antibacterial cleaning fluids and sprays
- The keyboard can be disabled during cleaning with the push of a button
- Automatic reminder for required cleaning (time-based, adjustable)*
- Three-level button illumination
- Medium button movement for an optimal typing experience
- Wear-resistant, removable surface made using a smooth, silicone-based material
- Multiple keyboard layouts for different languages

We Listen. We Adapt. We Deliver.

At Reinsberg Group, our mission is to ignite innovation, fuel collaboration, and pursue excellence across our dynamic portfolio in our Group companies.

We are dedicated to revolutionizing healthcare solutions, championing sustainable growth, and delivering exceptional value to our stakeholders.

All this while upholding the highest standards of integrity and social responsibility.

Rein Medical has been a member of the Reinsberg Group since 2026.

Member of REINSBERG® GROUP

We Listen. We understand what our customers truly need: For over two decades, Rein Medical has been developing customized IT solutions for the healthcare sector. Through close collaboration with medical practices, clinics, and partners, we analyze requirements, identify opportunities, and lay the groundwork for sustainable innovation. Our goal is not merely to deliver solutions, but to create genuine added value in everyday clinical practice.

We Adapt. We are flexible, personalized, and forward-thinking. As a systems service provider, we support projects holistically: from planning and development through to implementation. In Mönchengladbach, we develop high-quality all-in-one computer solutions, displays, IT wall modules for operating rooms, and intelligent video, audio, and image management systems. Our strength lies in the combination of in-house hardware and software development and our ability to tailor solutions precisely to our customers' needs.

We Deliver. We deliver—reliably, precisely, and to the highest quality standards. For us, “Made in Germany” is not just a promise, but a standard we live by. With our ISO 13485 certification, we ensure consistent quality management and continuous improvement. Together with a strong network of international partners, we ensure that our solutions guarantee smooth, efficient, and safe procedures in the operating room—today and in the future.

Headquarter in Germany
Rein Medical GmbH
Monforts Quartier 23
41238 Mönchengladbach
Germany
Tel.: +49 2161 6984-0
Fax: +49 2161 6984-259
E-mail: info@reinmedical.com

Switzerland
Rein Medical AG
Büfelderstrasse 1
8370 Sirmach TG
Switzerland
Tel.: +41 71 929-55-99
Fax: +41 71 929-55-90
E-mail: info.ch@reinmedical.com

Spain
Rein Medical Systems S.A.
C/ Téllez, 30
1a Planta, Oficina 1 Puerta 2
28007 Madrid, Spain
Tel.: +34 91 530-88-24
Fax: +34 91 574-32-93
E-mail: info.es@reinmedical.com